


➤ Hei-VAP Value / Hei-VAP Value Digital

Superior ease of use due to direct control

The Hei-VAP devices with hand lift are optimally suited for all standard applications. Controlling parameters is done so with the simple turn of a knob

- The convenience of large dial controls for adjustment of rotation speed and bath temperature
- **Hei-VAP Value Digital** models come with a digital display for set and actual heating bath temperature

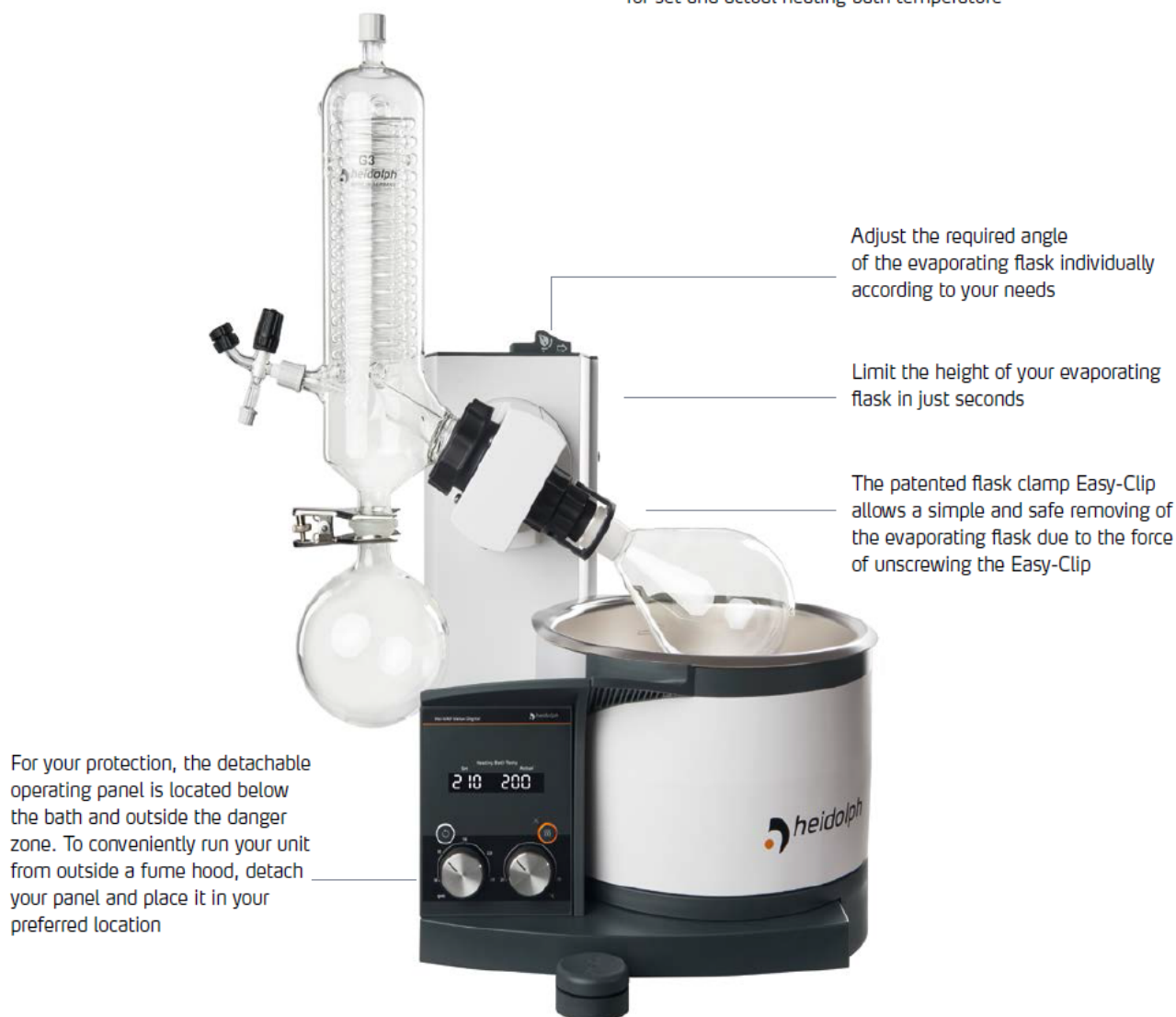

Hei-VAP Value Digital HL/G3
P/N 560-01302-00

Model	Hei-VAP Value hand lift (HL)		Hei-VAP Value Digital hand lift (HL)	
	standard	coated	standard	coated
G1 Diagonal	560-01100-00	560-01110-00	560-01102-00	560-01112-00
G3 Vertical	560-01300-00	560-01310-00	560-01302-00	560-01312-00
G3 XL Vertical	560-01305-00	560-01315-00	560-01307-00	560-01317-00
G5 Dry Ice Condenser	560-01500-00	560-01510-00	560-01502-00	560-01512-00
G6 Vertical for Reflux	560-01600-00	560-01610-00	560-01602-00	560-01612-00

➤ Glassware

- Glassware set includes one 1,000-ml evaporating flask and one 1,000-ml receiving flask
- All glassware sets are also available with transparent plastic coating for added safety
- Evaporating flasks and vapor tubes which come with a standard NS 29/32 joint size are also available with NS 24/29 ground joints
- Also available upon request: glassware sets for Soxhlet extraction and descending condenser system
- All glassware sets feature GL 10 thread; the inlet tube can be connected with a threaded fitting
- Glassware sets G3 and G6 allow for the automated distillation program AUTO_{accurate} on Hei-VAP Precision models

XL Condenser


For maximum evaporation performance

The new XL condenser with industry-leading condensing surface of 2,200 cm² achieve up to 40 % higher evaporation rates. Available in G3 XL and G9* XL versions




Glassware set G1

- Diagonal condenser for all standard distillations; the option with the **most affordable price**

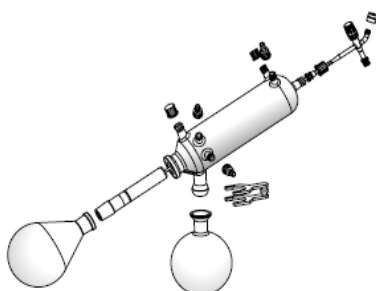


Glassware set G3

- Vertical condenser for all standard distillations; the option with the smallest space requirements, also available as XL version

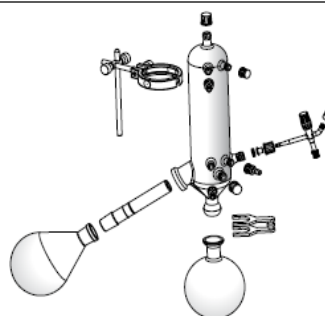


Glassware set G5

- **Dry Ice** condenser for low-boiling solvents

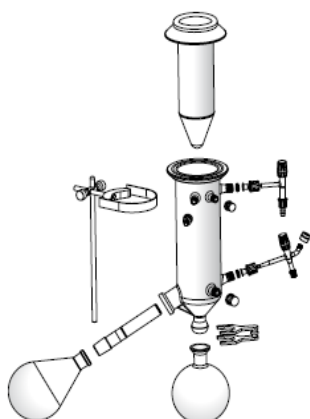

Glassware set G6

- Vertical condenser; centerpiece with valve for **reflux distillation**

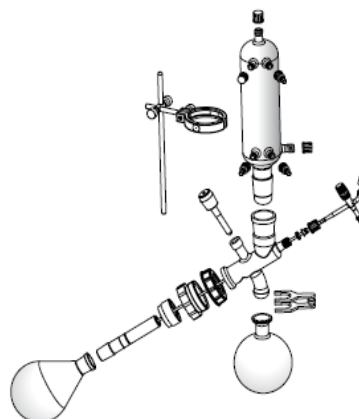

* For use with distimatic bench-top