

SF-110B Series

Impact Printer



A Trusted Printer that Supports Industry Compliance

With a built-in real time clock and user management, the OHAUS SF-110B series Impact Printer ensures data integrity, and meets GLP/GMP requirements. An easy-to-read OLED display and advanced functions make the SF-110B ideal for use with all OHAUS balances and scales.

Unique Features:

- A login password for each user fully prevents unauthorized changes to the printer settings such as date/time.
- The SF-110B provides users with full functionality from the most basic printing needs to the advanced applications, including statistics, totalization, formulation, averaging and interval printing.
- Additional information like Sample name/ID can be printed before each weighing that complements the report and improves data traceability.

SF-110B Series Impact Printer

Model	SF-110B	
Printing method	8 print solenoids impact dot matrix.	
Printing unit	Dot matrix printer 5x7, 40-character line length	
Printing speed	1.7 lines per second	
Lifetime	1 million lines	
Operation modes	Direct, Weighing, Multi Weighing, Statistic, Sum, Mean, Formulation, Loss Weight	
Paper roll (W x Dia.)	Standard paper, 57.5 mm x Ø 50 mm max.	
Ink ribbon cartridge	EPSON ERC-22 or ERC-09B	
Display	128 x 64 dots, OLED display	
Power Requirements	AC adapter Input: 100 – 240 VAC, 50 - 60 Hz AC adapter output: 12VDC, 2.0A	
Interface	RS232C	
Ambient temperature	0°C – 40°C	
Printer Dimensions (W x D x H)	115.00 mm x 191.27 mm x 88.32 mm	
Packaging dimensions (W x D x H)	177 mm x 235 mm x 169 mm	
Net weight	580g (paper roll included)	
Shipping Weight	1.5kg	

Compliance

- Product Safety: IEC/EN 62368-1
- Electromagnetic Compatibility: EN 55032; EN 55035, EN 61000-3-2, EN 61000-3-3
- Compliance Marks: CE, RCM, UKCA
- Other: WEEE, RoHS

Other Standard Features


High-speed dot matrix printing, 40 characters per line, auto tare operation directly from the printer, multi-user access control, editable header and trailer. Fully compatible with OHAUS weighing instruments.

Multi-Language printing in DIRECT mode: English, Spanish, German, Italian, French, Chinese, Korean, Japanese, Polish, Portuguese, Hungarian, Turkish and Czech.

Accessories

Paper Roll	12120799
Ink Ribbon, Black	30529322

Outline Dimensions

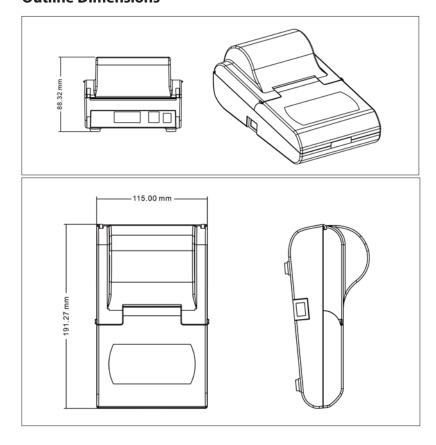

OHAUS Europe GmbH

Heuwinkelstrasse 3, 8606 Nänikon, Switzerland

e-mail: ssc@ohaus.com Tel: 0041 22 567 53 19 e-mail: tsc@ohaus.com Tel: 0041 22 567 53 20

www.ohaus.com

The management system governing the manufacture of this product is ISO 9001:2015 certified.



80776695_A 20220530 © Copyright OHAUS Corporation