Model CBF 260 | CO $_{2}$ incubators with hot air sterilization and humidity regulation

The BINDER CBF is the high-end CO₂ incubator with humidity control that is suitable for all sensitive incubation applications in GMP environments. It ensures consistent results including in cell-based assays with multiwell plates. This BINDER incubator is even suitable for complex cultivation experiments or individual environments with hypoxic conditions thanks to its comprehensive assortment of options and accessories.

BENEFITS

- Safe cell cultures receive maximum protection against contamination
- Reliable constant well-being conditions for the cells
- Smart simple routine cleaning and convenient operation
- Economical efficient operation without consumables





Model 170

MAIN FEATURES

- Contamination-safe humidification system without water pan
- Tried-and-tested anti-contamination concept with 180 °C hot air sterilization
- Hot-air sterilizable CO₂ sensor with single-beam infrared technology
- Seamless inner chamber made of stainless steel with flanges as shelf support system
- Fanless interior with Venturi CO₂ gas mixing nozzle
- Intuitive touchscreen controller
- Internal data logger, measured values can be read out in open format via USB
- Troubleshooting system with visual and acoustic alarms
- Interfaces: Ethernet, USB, zero-voltage alarm contact
- Fail-safe CO2 system for protecting the pH of cell cultures

Model 170

OPTIONAL EQUIPMENT

- O_2 control with O_2 control ranges: 0.2 to 20 vol.% O_2 or 10 to 95 vol.% O_2
- Multi-sectional glass doors with special shelves for dividing up the space to accommodate different experiments
- Stacking adapter for secure, space-saving installation
- Gas tank changer for connecting two gas tanks to one or two incubators

ORDERING INFORMATION

Interior volume [L]	Voltage	Option model	Version	ArtNo.
	200230 V 1~ ph 50/60 Hz	Standard	CBF260-230V	9640-0021
		with O_2 control 0.2–20%	CBF260-230V-0	9640-0022
267		with O_2 control 0.2–20% and cell therapy compartmentalization	CBF260-230V-OZ	9640-0025
		Standard	CBF260UL-120V	9640-0023
	100120 V 1~ ph 50/60 Hz	with O_2 control 0.2–20%	CBF260UL-120V-0	9640-0024
		withO ₂ control 0.2–20% and cell therapy compartmentalization	CBF260UL-120V-OZ	9640-0026



TECHNICAL DATA

Description	CBF260-230V ¹	CBF260-230V-01	CBF260-230V-0Z ¹	CBF260UL-120V ¹	CBF260UL-120V- 01	CBF260UL-120 OZ ¹
Article Number	9640-0021	9640-0022	9640-0025	9640-0023	9640-0024	9640-0026
Performance Data Temperature						
Temperature range +6 °C above ambient temperature to [°C]	60	60	60	60	60	60
Temperature variation at 37 °C [± K]	0.4	0.4	0.4	0.4	0.4	0.4
Temperature fluctuation at 37 °C [± K]	0.1	0.1	0.1	0.1	0.1	0.1
Recovery time after 30 seconds door open at 37 °C [min]	6	6	6	6	6	6
Performance Data Climate						
lumidity range	5095	5095	5095	5095	5095	5095
CO ₂ -Data						
CO ₂ range [Vol% CO ₂]	020	020	020	020	020	020
CO ₂ measurement	IR	IR	IR	IR	IR	IR
CO_2 Recovery time after 30 seconds door open at 5 Vol% CO_2 [min]	5	5	5	5	5	5
D₂-Data						
Standard O₂ control ranges: O₂ range [Vol% O₂]		0.220	0.220		0.220	0.220
O_2 control ranges with option: O_2 range [Vol $\%$ O_2]		1095	1095		1095	1095
D₂ Recovery time after 30 seconds door open at 5 Vol% O₂ [min]		15	15		15	15
Electrical data						
Rated Voltage [V]	200230	200230	200230	100120	100120	100120
Power frequency [Hz]	50/60	50/60	50/60	50/60	50/60	50/60
Nominal power [kW]	1.5	1.5	1.5	1.5	1.5	1.5
Jnit fuse [A]	10	10	10	16	16	16
Phase (Nominal voltage)	1~	1~	1~	1~	1~	1~
Measures						
nterior volume [L]	267	267	267	267	267	267
Net weight of the unit (empty) [kg]	125	125	125	125	125	125
Permitted load [kg]	40	40	40	40	40	40
Load per rack [kg]	10	10	10	10	10	10
Nall clearance back [mm]	100	100	100	100	100	100
Wall clearance sidewise [mm]	50	50	50	50	50	50
nternal Dimensions						
nterior width [mm]	620	620	620	620	620	620
nterior height [mm]	750	750	750	750	750	750
nterior depth [mm]	575	575	575	575	575	575
	575	575	575	575	575	575
Doors	6	575 6	575 6	575 6	575 6	575 6
Doors nner doors						
Interior depth [mm] Doors Inner doors Unit doors Housing dimensions not incl. fittings and	6	6	6	6	6	6
Doors nner doors Jnit doors Housing dimensions not incl. fittings and	6	6	6	6	6	6
Doors Inner doors Jnit doors	6 1 I connections	6 1	6	6	6	6


1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

Data Sheet Model CBF 260

Description	CBF260-230V ¹	CBF260-230V-01	CBF260-230V-0Z1	CBF260UL-120V ¹	CBF260UL-120V- 01	CBF260UL-120V- OZ ¹
Article Number	9640-0021	9640-0022	9640-0025	9640-0023	9640-0024	9640-0026
Environment-specific data						
Energy consumption at 37 °C [Wh/h]	65	65	65	65	65	65
Sound-pressure level [dB(A)]	41	41	41	41	41	41
Fixtures						
Number of shelves (std./max.)	3/8	3/8	3/8	3/8	3/8	3/8

1 All technical data is specified for unloaded units with standard equipment at an ambient temperature of +22 °C ±3 °C and a power supply voltage fluctuation of ±10%. The temperature data is determined to BINDER factory standard following DIN 12880, observing the recommended wall clearances of 10 % of the height, width, and depth of the inner chamber. All indications are average values, typical for units produced in series. We reserve the right to change technical specifications at any time.

DIMENSIONS incl. fittings and connections [mm]

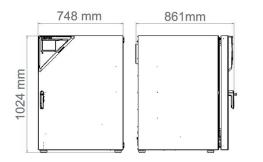

OPTIONS

Designation	Description	*	ArtNo.
Access port with silicone plug	back		
	30 mm	01	8612-0025
	left		
	30 mm	01	8612-0026
	right		
	30 mm	01	8612-0027
Analog output 4-20 mA	for temperature and CO ₂ values (outputs not adjustable)	02	8612-0022
Calibration certificate O ₂	for versions with O_2 control; a measuring point in center of chamber at a specified O_2 percentage value	-	8012-0229
	for temperature, measurement in center of chamber at specified temperature	-	8012-1132
Calibration certificate,	temperature measurement incl. certificate and 27 measuring points at specified temperature	-	8012-1592
temperature	temperature measurement incl. certificate, 15- 18 measuring points at specified temperature	-	8012-1571
	temperature measurement incl. certificate, 9 measuring points at specified temperature	-	8012-1550
Calibration certificate, temperature and CO ₂	for temperature and CO_2 , temperature measurement in center of chamber, CO_2 measurement performed using test gas at 37 °C and 5 % CO_2	-	8012-1235
CELLROLL Set	modular, expandable roller bottle system for cell cultivation, 6 roller bottles	-	8612-0024
Door hinged on the left	Unit door and standard glass inner door with door hinged on the left	-	8612-0034
Electric access port	8-pin, for low voltage with LEMO socket (coverable) and LEMO plug (max. 24 V – 2 A)	-	8612-0033
Inner door, divided	Cell therapy compartmentalization, consisting of 8 inner doors and four shelf levels with one partition wall each	-	8612-0031
O ₂ control range	alternative to the standard range, for controlling O_2 concentrations between 10 and 95 vol. %, only in conjunction with option models with O_2 control	-	8612-0028
Quick sample access	vewing window (13 x 25 cm) in the inner door for quick access	-	8612-0032

* Notes > See last page



ACCESSORIES

Designation	Description	*	ArtNo.
APT-COM™ 4 BASIC-	for simple logging and documentation requirements with up to 5 networked units.		
dition	version 4, BASIC edition	-	9053-0039
APT-COM™ 4 GLP-	for working under GLP-compliant conditions. Measured values are documented in a tamper-proof way in line with the requirements of FDA Regulation 21 CFR 11.		
Edition	version 4, GLP edition	-	9053-0042
APT-COM™ 4	convenient unit and user management built on the BASIC edition. Suitable for networking up to 100 units.		
PROFESSIONAL-Edition	version 4, PROFESSIONAL edition	-	9053-0040
Base	the base equipped with casters is used for the secure positioning and leveling of a BINDER $\rm CO_2$ incubator	-	9051-0044
Convenient water supply	Set consisting of 10 liter water container with magnetic holder, hose adapter with practical quick coupling and water stop function, flexible positioning, enables water to be supplied to multiple CB units simultaneously via a water container, not pre-filled.		
	for versions with active humidification	-	8009-1128
	for CO_2 , consisting of a gas tank pressure regulator with connection parts and 5-meter hose	-	8012-0014
Gas cylinder connection	for N_2 , consisting of a gas tank pressure regulator with connection parts and 5-meter hose	-	8012-0016
	for O ₂ , consisting of a gas tank pressure regulator with connection parts and 5-meter hose	_	8012-0015
pH-neutral detergent	concentrated, for gentle remove of residual contaminants; 1 kg	-	8012-2250
	IQ/OQ documents – supporting documents for validation performed by customers, consisting of: IQ/OQ checklists incl. calibration guide and comprehensive unit documentation; parameters: temperature, CO ₂ , O ₂ , pressure, depending on unit		
	Digital in PDF format	-	7057-0001
Dualification de sumante	Hard copy inside folder	-	7007-000
Qualification documents	$IQ/OQ/PQ$ documents – supporting documents for validation performed by customers, according to customer requirements, PQ section added to qualification folder IQ/OQ; parameters: temperature, CO_2 , O_2 – or pressure, depending on unit		
	Digital in PDF format	-	7057-0005
	Hard copy inside folder	-	7007-000
Rubber pads	set anti-slip feet	-	8012-2030
Shaker "Professional"	Shaker and controller equipped with LEMO connectors, exclusively for CB with electric access port, 100 – 240 V, 50/60 Hz	-	8012-1772
Shaker "Universal"	Shaker and controller connected via ribbon cable, universally applicable, 100 – 240 V, 50/60 Hz	-	8012-1771
	Stainless steel	-	8012-2183
Shelf, perforated	for divided inner door		
	Stainless steel, divided	-	8012-1996
Stacking adapter	for the thermally isolated stacking of two BINDER CO $_{\scriptscriptstyle 2}$ incubators	-	9051-0039
Stacking frame	stable, vibration-free stacking frame on casters with stop brake for safe stacking of two CB series units	-	9051-0041
	Set consisting of 10 bottles, each with 1,000 ml of sterile humidification water, including peristaltic tubing, connectors, and fine gas filters		
	for versions with active humidification	-	8009-1104
Water supply	Set consisting of 12 empty, refillable tubular bags including peristaltic tubing, connectors, and fine gas filters		
	for versions with active humidification	-	8012-2255
	Set consisting of 3 empty, refillable tubular bags including peristaltic tubing, connectors, and fine gas filters		
	for versions with active humidification	-	8012-2256
WLAN kit	The kit contains a client bridge that establishes a wireless connection between BINDER units and APT-COM4, LIMS, or customer-specific software via an Ethernet interface. It provides an alternative solution in situations where units need to be positioned in locations without an on-site Ethernet connection. In secured networks, installation and configuration must be performed by the customer's IT service.	-	8012-2262

* Notes > See last page

SERVICES

Designation	Description	*	ArtNo.
Maintenance contracts			
nstallation services			
Maintenance services			
Calibration services			
Validation services			
Warranty service			
1-year warranty extension	The warranty is extended by 1 year from the delivery date, wear parts are excluded	-	DL50-0020
BRONZE 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts	05	DL20-0710
Colibration cartificate temporature	Calibration of one (1) test temperature specified by the user in center of chamber, including certificate	03, 04, 05	DL30-0101
Calibration certificate, temperature	Extension of calibration of one (1) additional test temperature specified by the user in the center of the usable space, including certificate	03, 04, 05	DL30-010:
CO ₂ calibration	CO_2 calibration at a specified value, measurements are taken with analyzed test gas at 5%, including certificate	03, 04, 05	DL30-040
Execution of IQ/OQ	Execution of IQ/OQ in accordance with qualification folder	05	DL41-0210
Execution of IQ/OQ/PQ	Execution of IQ/OQ/PQ in accordance with qualification folder	05	DL44-050
GOLD 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, replacement of wear parts, calibration of one temperature and CO_2 value, including certificate	05	DL20-094
Maintenance	One-off maintenance service in accordance with maintenance schedule. Visual inspection of mechanical and electrical components, testing of all key functions. Calibration of a test temperature specified by the user in center of usable space without certificate	05	DL20-040
O₂ calibration	$O_{\rm 2}$ calibration with 1 measuring point in center of chamber at a specified $O_{\rm 2}$ percentage value, including certificate	03, 04, 05	DL30-040
SILVER 3-year maintenance contract	Maintenance service as contractually agreed, visual inspection of mechanical and electrical components, check of control response, 20% discount on spare parts, testing of all key functions, calibration of one test temperature specified by the user in the center of the usable space, without certificate	05	DL20-082
Temperature measurement, 18 measuring points	Temperature measurement with 18 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-0118
Temperature measurement, 27 measuring points	Temperature measurement with 27 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-012
Temperature measurement, 9 measuring points	Temperature measurement with 9 measuring points with a set value specified by the user, including certificate	03, 04, 05	DL30-010
Jnit commissioning	Connect the unit to the customer-side connections (electricity, water, wastewater, gas), basic functions check, brief operating instructions. (excl.: unpacking, setup, controller instructions, programming, installation work)	05	DL10-010
Unit instructions	Instruction regarding operating principle and basic functions of the unit, operation of the control electronics including programming	05	DL10-050

* Notes › See last page

NOTES

- Condensation may occur in the area around the access port. Access ports may be placed in custom locations for an additional charge. UL mark is not granted when this option is used. 01
- 02
- Sensor calibration is performed in an accredited calibration laboratory. 03
- 04 Calibration is performed according to the BINDER factory standard.
- 05 Quoted prices do not include travel costs. Please refer to the chapter on BINDER Service for travel costs for your region. Quoted prices for services performed in Switzerland do not include a country-specific added fee (available on request).

BINDER GmbH

Tuttlingen, Germany TEL +49 7462 2005 0 FAX +49 7462 2005 100 info@binder-world.com www.binder-world.com

BINDER Asia Pacific (Hong Kong) Ltd.

Kowloon, Hong Kong, P.R. China TEL +852 39070500 FAX +852 39070507 asia@binder-world.com www.binder-world.com

BINDER Environmental Testing

Equipment (Shanghai) Co., Ltd. Shanghai, P.R. China TEL +86 21 685 808 25 FAX +86 21 685 808 29 china@binder-world.com www.binder-world.com

BINDER Inc.

Bohemia, NY, USA TEL +1 631 224 4340 FAX +1 631 224 4354 usa@binder-world.com www.binder-world.us