

## Product data sheet

# VALEGRO 350 Recirculating Cooler


## Modern, User-friendly Chillers for a Wide Range of Cooling Tasks

The air-cooled VALEGRO recirculating coolers are environmentally friendly all-rounders for a wide range of cooling and heating tasks from -20 to +40 °C. The pump can be user-friendly adjusted in different pump stages on a percentage basis.


Whether in the laboratory, research, or production, the VALEGRO 350 recirculating coolers offer a constant cooling capacity of 0.35 kW at 20 °C, supporting stable processes around the clock. These modern devices offer precise temperature control for a wide range of applications and, thanks to their compact design, fit seamlessly into any working environment. The refrigeration circuit uses an eco-friendly, natural refrigerant. When designing the VALEGRO recirculating coolers, the focus was on maximum user-friendliness. Many helpful features such as a timer function, a removable venting grid that can be removed without tools, and an illuminated and simple to clean level indicator make daily work easier. The inclined OLED touch display is easy to read and ergonomic to use. The intuitive menu navigation allows quick access to all functions. Warning notices and system messages are displayed in clear, understandable language, ensuring maximum safety when using the device. In addition to RS232 and USB-C interfaces, VALEGRO recirculating coolers are also available with Ethernet (optional). All interfaces are located directly below the display and are easy to use from the front.




Pump capacity

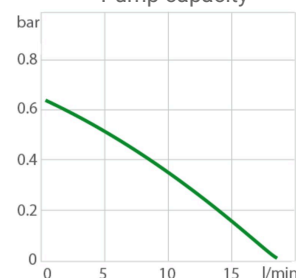

Heat-up time Medium: Thermal G

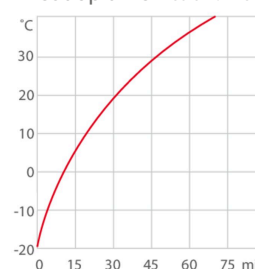

Cool-down time Medium: Ethanol

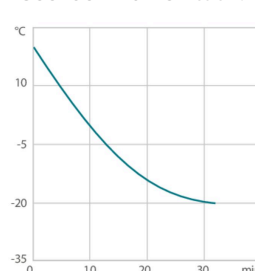

### Product Features

- ✓ Precise temperature stability of  $\pm 0.3$  °C
- ✓ Removable venting grid without tools
- ✓ Integrated castors for good handling on and under the lab bench
- ✓ Simple timer function
- ✓ Autostart function
- ✓ Data recording
- ✓ Intuitive menu navigation in full text German/English
- ✓ User-friendly operation thanks to comprehensive help texts

## Available options

### Interfaces

#### Activation Ethernet interface

The Ethernet interface integrated in the VALEGRO recirculating cooler is activated at the factory and enables the use of network-based communication and control functions.



## Performance values

230V/50Hz (Schuko Plug - CEE 7/4 Plug Type F)

Viscosity max. cSt	50
Pump capacity flow pressure l/min	16
Pump capacity pressure bar	0.6
Power consumption A	3
Voltage Tolerance %	230V   ±5

Order number including voltage version (230V/50Hz) 9610035.30003787

#### Cooling capacity 1 (Ethanol)

°C	20	0	-10	-20
kW <sup>1</sup>	0.35	0.25	0.18	0.07

Note about natural refrigerants:

Temperature control units using natural refrigerants are often subject to regulatory requirements regarding the installation site, operation, transport or disposal of the units. If you have any questions, we will be happy to advise you.

#### Refrigerant stage 1

Refrigerant	R290
Filling weight g	35
Global Warming Potential for R290	3
Carbon dioxide equivalent t	0.000105

<sup>1</sup> Performance specifications measured in accordance with DIN 12876. Cooling capacities up to 20 °C measured with ethanol; over 20 °C with thermal oil unless otherwise specified. Performance specifications apply at an ambient temperature of 20 °C. Performance values may differ with other bath fluids.

## Technical data

### Available voltage versions

Order No. 9610035

### Available voltage versions:

9610035.30003785 100V 50/60Hz (Nema N5-15 Plug) (R290)

9610035.30003786 115V 60Hz (Nema N5-15 Plug) (R290)

9610035.30003787 200-230V 50/60Hz (Schuko Plug - CEE 7/4 Plug Type F) (R290)

9610035.30003788 230V 50Hz (UK Plug Type BS1363A) (R290)

9610035.30003789 230V 50/60Hz (CH Plug Type SEV 1011) (R290)

9610035.30003790 200-230V 50/60Hz (CN Plug) (R290)

### Cooling

Cooling of compressor

1-stage Air

### Bath

Bath tank

Stainless steel

### Electronics

Interfaces

Ethernet optional, RS232, USB

Temperature control

PID1

Temperature display

OLED

Temperature setting

Touchscreen

Electronic Timer h:min

99 ... 59

### Dimensions and volumes

Weight kg

26

Barbed fittings inner diameter mm

8/12 mm

Total dimensions cm (W x L x H)

23.5 x 45 x 45

Body dimensions cm (W x L x H)

23 x 40 x 44.5

Filling volume l

1.5 ... 3

Pump connections

M16x1 male

### Temperature values

Return flow temperature max. °C

+85

Working temperature range °C

-20 ... +40

Temperature stability °C

±0.3

Ambient temperature °C

+5 ... +40

Temperature display resolution °C

0.1

### Other

Classification

Classification I (NFL)

Standards

CE, RoHs

Alarm notification

Acoustic, Optical

IP Code

IP 21

Pump function

Pressure Pump

Pump type

Centrifugal Pump

User Interface Language

English, German

Sound pressure level dbA

52

#### Included in delivery

2 Barbed fittings for tubing 8 and 12 mm ID. (Pump connections  
M16x1 male)

Refer to [www.julabo.com](https://www.julabo.com) for more information regarding the entire JULABO product portfolio. Technical changes without prior notification.  
Images may deviate from the original. | Datasheet No.EN9610035/251218