

PRODUCT SPECIFICATION - TT300E50-5K

Environmental test chamber TESTA TT300E50-5K




TESTA TT – Precise and reliable environmental test chambers.

The TESTA CT and TESTA TT environmental test chambers expand the Memmert portfolio with high-performance solutions for environmental simulation.

With a wide setting temperature range of -52 to +180 °C and configurable temperature change rates of up to 10 K/min, the chambers offer an ideal solution for temperature and climatic testing, simulations of environmental conditions, ageing tests and stress tests.

Features

- Configurable temperature Change Rates: up to 10K/min
- State-of-the-Art Connectivity
- Flexible Interior Design
- Premium Quality

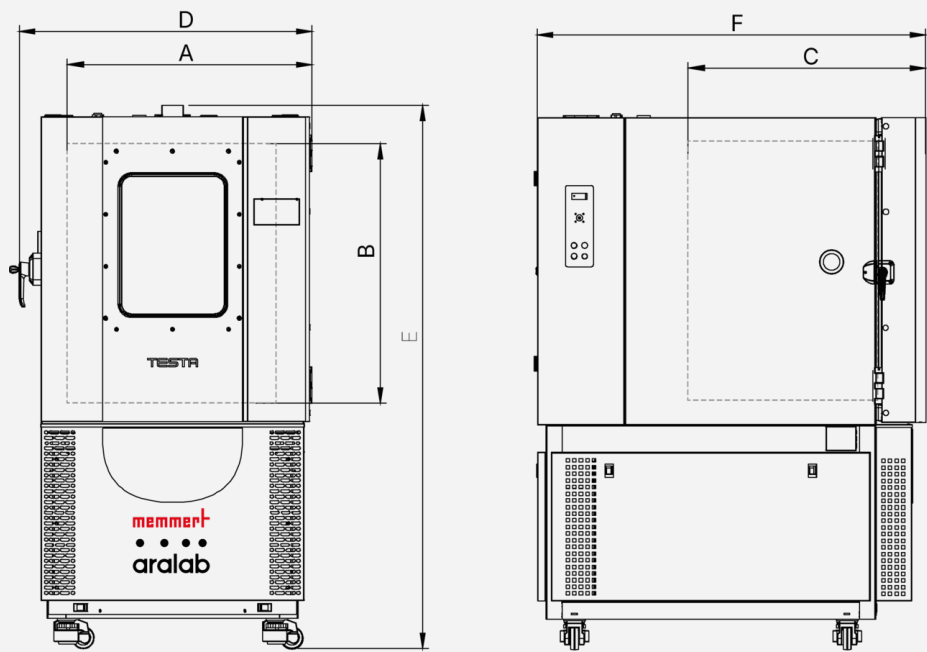

Temperature

Setting temperature range	°C	-52 to +180
Setting accuracy	°C	0.1
Temporal temperature deviation		± 0.1 to ± 0.3
Temperature control system		By stainless steel electric heaters located in the air treatment tunnel

Features

Control unit		Touch screen controller
Adjustable parameters		Temperature
Adjustable device functions		• Climate control regulation • Real-time display of process variables • Program execution • Alarm definition • and much more
Interfaces		Ethernet LAN

Safety

Parameter monitoring		1 Pt100 sensor DIN class A 1 Pt100 sensor DIN class A, movable sensors for flexible placing inside chamber
Temperature monitoring		Safety thermostat with High / Low temperature configuration, with automatic stop of all thermic systems
Alarm		High / Low Temperature alarms

Heating concept

Cooling system		Air cooled hermetic compressor group (low noise and high efficiency) with enforced ventilation and without CFC's
----------------	--	--

Standard equipment

Door		Lockable door, double silicone door gaskets and anti-condensation heating frames
Standard option		- Entryport Ø 80 mm, left - Door locking, mechanical
Type standard insert		Stainless steel grid, electropolished
Number of slide-in units		2
Other standard accessories		- Water conductivity measuring device - Documentation package

Stainless steel interior

Volume	l	272
Width	mm	634
Height	mm	715
Depth	mm	600
Max. number of grids/shelves	number	10
Max. loading of chamber	kg	120
Max. loading of grids/shelves	kg	25

Textured stainless steel casing

Width	mm	1009
Height	mm	1840
Depth	mm	1340
Installation		on height-adjustable and lockable castors

Electrical data

Electrical load; 400 V or 3x230 V	kW	11
-----------------------------------	----	----

Ambient conditions

Distance to neighbouring appliances / wall clearances	back	mm	1000
	lateral	mm	500

Packing/shipping data

Transport information			The appliances must be transported upright
Customs tariff number			8419 8998
Country of origin			PT
Packaging	Width	mm	1180
	Height	mm	2100
	Depth	mm	1575
Weight	net	kg	513
	gross	kg	704

Standard units are safety-approved and bear the test marks:

