



EduBlue

HIGHLIGHTS

- Designed for education
- Robust and sturdy, with secured eyepieces
- WF 10x/20 mm eyepieces
- Dual or triple magnifications
- Binocular and trinocular models available
- Provided with APL (Anti-bacterial Protection Layer)
- Quick focus scale
- 1 W LED illumination
- With rechargeable batteries for cordless use
- Digital models with 4.0 MP WiFi camera available
- Ergonomic carrying grip
- 5 years warranty



TECHNICAL SPECIFICATIONS

EYEPIECE

Pair of secured WF10x/20 mm eyepieces supplied with eyecups

HEAD

- Binocular head with 45° inclined tubes
- Trinocular head with 45° inclined tubes and 23,2mm tube. The standard trinocular tube accepts all standard 23.2 mm C-mount objectives
- Diopter adjustment of ± 5 on one side
- Interpupillary distance adjustable between 55 and 75 mm
- Digital head is supplied with a 4.0 MP USB-2, 1/1.8" CMOS WiFi camera with a maximum of resolution of 2688 x 1512 pixels

DUAL MAGNIFICATION OBJECTIVES

- Revolving nosepiece with dual 1x/3x and 2x/4x, which can provide standard magnifications of 10x and 30x or 20x and 40x
- Working distance 60 mm
- Field of view of 20/6.7 mm or 10/5 mm
- Magnifications can be altered using optional eyepieces

All optics are anti-fungus treated and anti-reflection coated for maximum light throughput

TRIPLE MAGNIFICATION OBJECTIVES

- Revolving nosepiece with triple 1x/2x/3x and 1x/2x/4x, which can provide standard magnifications of 10x, 20x and 30x or 10x, 20x and 40x
- Working distance 60 mm
- Field of views of 20/10/6.7 mm or 20/10/5mm
- Magnifications can be altered using optional eyepieces

STAND

- The rack & pinion and pillar stands of the EduBlue are equipped with ergonomic designed flat bases, complete with 2 object clamps and $\varnothing$ 60 mm transparent and black/white stage plate
- The coarse adjustment is equipped with tension control
- The stands are alloy metal casted with hardened coating

ED.1302-P / ED.1402-P

225 (h) x 175 (w) x 220 mm (d) | 2.4 kg

ED.1302-S / ED.1402-S

260 (h) x 175 (w) x 220 mm (d) | 2.1 kg

ED.1303-S / ED.1403-S

260 (h) x 175 (w) x 220 mm (d) | 2.2 kg

ED.1405-S / ED.1805-S

260 (h) x 175 (w) x 220 mm (d) | 2.3 kg

ED.1302-EVO / ED.1402-EVO

260 (h) x 175 (w) x 220 mm (d) | 2.5 kg

- The rack & pinion stands of the EduBlue are equipped with graduations and focus point indicator for quick focussing
- EVO models are equipped with a build-in moveable object stage

DUAL ILLUMINATION SYSTEM


Easily switch between:

- 1 W white LED: 6000–6500 K, ideal for standard observations with crisp, cool illumination
- 1 W colored LED: 3000–3200 K, enhances contrast for specific applications, such as observing transparent objects
- 60 mm working distance
- Evo models feature an additional flexible arm with 1 W LED illumination

All EduBlue microscopes are supplied with rechargeable AA batteries for cordless operation and an external 5 V DC power adapter

PACKAGE CONTENT

- Supplied with an external 100-240 V mains adapter/charger and three rechargeable batteries for corded and cordless use, dust cover, $\varnothing$ 60 mm transparent and black/white stage plates, eyecups and user manual
- All packed in a polystyrene box

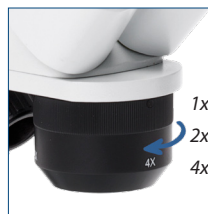

Triple magnification



Quick focus scale



Moveable object stage

DIGITAL MODELS

CAMERA

- Digital head is supplied with a 4.0 MP USB-2 1/1.8" CMOS WiFi camera
- A maximum resolution of 2688 x 1512, 30 fps (Frame rates have been determined at optimal conditions)
- 12 bits grayscale conversion and 24 bits color rendering
- The camera streams H.264 encoded high resolution images simultaneously to several WiFi enabled smartphones, tablets or computers with iOS, Android and Windows operating systems
- It can also operate via USB 2 and ethernet in combination with a PC or Mac computer
- Delivered with a micrometer 1mm/100 calibration slide
- Warranty for the camera is 2 years



APL (ANTIMICROBIAL PROTECTION LAYER)

Using the same technology already used extensively throughout the healthcare environment in wall coatings, patient lighting, air pressure stabilization and suction liners, APL prevents the growth of unwanted microbes that can cause degradation, discoloration, staining or odors of the Euromex microscopes

SOFTWARE

- ImageFocus Alpha software can be used for Windows based devices and WiFi connected PC and is supplied on CD
- Apps for Android for this camera are supplied on cd or can be downloaded from this page under Downloads (free of charge) or from the Apple Store (free of charge)
- Delivered with ImageFocus Alpha for capturing of images and videos
- This software also allows a full range of analysis like measurements on still and live images and annotations on captured images
- Save images in .jpg, .tif, .bmp or .dicom formats, save videos in .avi format
- Images can be annotated and measurements can be performed on live or captured images
- Compatible with Windows 10 and higher, all 32 and 64 bits configurations
- A Mac OS version is available
- Updates can be downloaded on our website www.euromex.com

MODELS	Bino	Trino	Digital head + WiFi	1x/3x objectives	2x/4x objectives	1x/2x/3x objectives	1x/2x/4x objectives	Pillar stand	Rack & pinion stand	Movable stage/ flexible light arm
ED.1302-P	•			•				•		
ED.1302-S	•			•					•	
ED.1303-S		•		•					•	
ED.1402-P	•				•			•		
ED.1402-S	•				•				•	
ED.1403-S		•			•				•	
ED.1405-SW			•		•				•	
ED.1502-S	•					•			•	
ED.1802-S	•						•		•	
ED.1805-SW			•				•		•	
ED.1302-EVO	•			•					•	•
ED.1402-EVO	•				•				•	•



● ED.1302-P



● ED.1405-S



● ED.1303-S



● ED.1402-EVO

ACCESSORIES AND SPARE PARTS

- ED.6005** Pair of HWF 5x/22 mm eyepieces
ED.6010 Pair of HWF 10x/20 mm eyepieces
ED.6015 Pair of HWF 15x/12 mm eyepieces
ED.6020 Pair of HWF 20x/10 mm eyepieces
ED.6110 HWF 10x/20 mm eyepiece with 10 mm/100 micrometer
ED.6099 Pair of eyecups
ED.9570 Pair of object clamps for stage
ED.9950 Stage plate frosted glass, opaque, Ø 60 mm
ED.9956 Stage plate black/white, Ø 60 mm
ED.9975 External 100-240 V mains adapter/ charger
ED.4300 Aluminium transport case for EduBlue microscopes
 Not suitable for ED.3000

- AE.9918** Nylon microscope bag (25 x 39 x 19 cm)
AE.9912 Euromex dust cover medium (35 x 39 cm)
SL.5504 LED replacement unit for EduBlue, incident illumination
SL.5505 LED replacement unit for EduBlue, transmitted illumination
AE.9918 Nylon microscope bag (25 x 39 x 19 cm)
AE.9912 Euromex dust cover medium (35 x 39 cm)
AE.1112 Micrometer 76 x 26 mm slide, 50 mm/50 divisions
PB.5245 Lens cleaning paper, 100 sheets per pack
PB.5274 Isopropyl alcohol 99% (200 ml)
PB.5275 Cleaning kit: lens cleaning fluid, lint free lens tissue, brush, air blower, cotton swabs



● ED.4300





● AE.9918



● AE.9912

DIMENSIONS

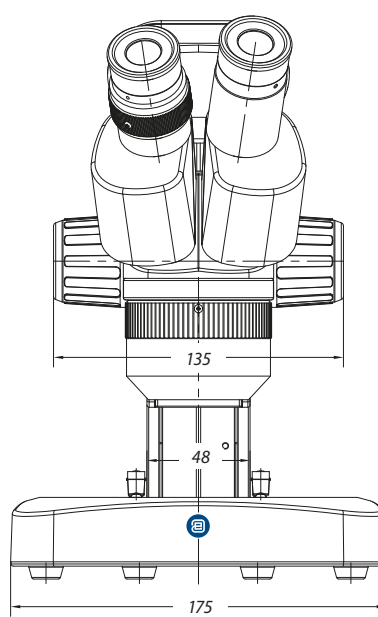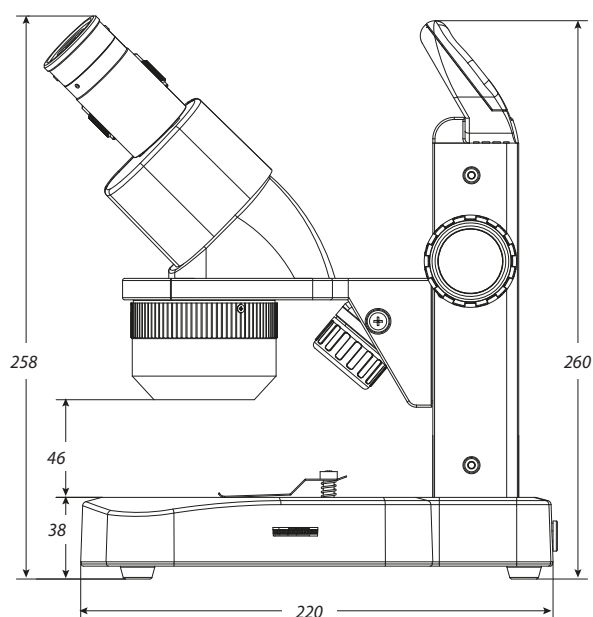

www.euromex.com