


Clean bench

WLC-V

Class 100 HEPA filter & vertical air flow

The laminar flow clean bench provides reliable protection for samples and work processes against particulate contamination. The air from above the cabinet vertically goes through a pre-filter and a class 100 HEPA filter with an efficiency up to 99.97% & 0.3 micron size.

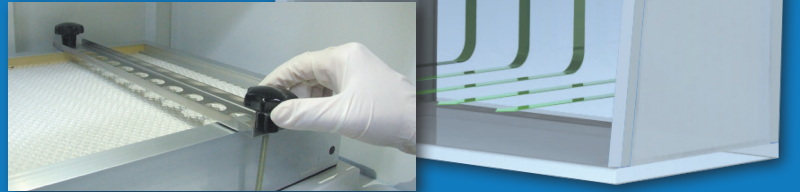

UV lamp

The built-in UV lamp contributes to a better sterilization. The UV lamp automatically turns off when the front door is opened.



Door

Easy to open tempered glass sliding door ensures safe operation and protection against UV radiation.



Operation of different types of devices within the clean bench due to two built-in power connections.



Connection of two different gas types possible.

Connections

Controller

Intuitive capacitive touch-button controller. 8-step control of air-flow speed. Timer function for UV lamp.

